# ONETOUCH<sup>®</sup> VerioIQ<sup>®</sup>

Blood Glucose Monitoring System

## Owner's Booklet

Instructions for use.

The system described herein is covered by one or more of the following U.S. patents: 6,179,979, 6,193,873, 6,284,125, 6,716,577, 6,749,887, 6,797,150, 6,863,801, 6,872,298, 7,045,046, 7,498,132, and 7,846,312. Use of the monitoring device included herein is protected under one or more of the following U.S. patents: 6,413,410, 7,749,371. Purchase of this device does not act to grant a use license under these patents. Such a license is granted only when the device is used with OneTouch® Verio™ Test Strip. No test strip supplier other than LifeScan is authorized to grant such a license. The accuracy of results generated with LifeScan meters using test strips manufactured by anyone other than LifeScan has not been evaluated by LifeScan.

Manufactured By:  
LifeScan Europe  
Division of Cilag GmbH  
International  
6300 Zug  
Switzerland

Meter made in China

If you have questions about the use of any OneTouch® product, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. - 10 p.m. Eastern Time) [www.OneTouch.com](http://www.OneTouch.com)

If you cannot reach Customer Service, contact your health care professional for advice.

© 2011 LifeScan, Inc.

Rev. date: 06/2014

Blood Glucose Monitoring System

## Owner's Booklet

# Owner's Booklet

## Thanks for choosing OneTouch®!

The OneTouch® Verio®IQ Blood Glucose Monitoring System is one of the latest product innovations from OneTouch®. Every OneTouch® Meter is designed to make it easy to test your blood glucose and help you manage your diabetes.

This Owner's Booklet offers a complete explanation of how to use your new meter and testing supplies, and reviews the do's and don'ts of testing your blood glucose level. Clear instructions and visuals will help you get comfortable using the OneTouch® Verio®IQ Meter, so you can expect to get reliable results every time you test. Please keep your Owner's Booklet in a safe place; you may want to refer to it in the future.

We hope OneTouch® products and services will continue to be a part of your life.

## Symbols

Cautions and Warnings. Refer to the Owner's Booklet and inserts that came with your system for safety-related information.

Direct current

Consult instructions for use

Battery full

Battery 30–85% full

Battery low

Battery very low

Before meal tag

After meal tag

Message icon

Control solution result

# Table of Contents

<b>1</b>	<b>Setting up your system</b>	<b>12</b>
<b>2</b>	<b>Taking a test</b>	
	- Testing your blood glucose	27
	- Testing with control solution	47
<b>3</b>	<b>Attaching tags to your blood glucose results</b>	<b>55</b>
<b>4</b>	<b>Reviewing past results, averages and high/low pattern messages</b>	<b>60</b>
<b>5</b>	<b>Care and maintenance</b>	<b>76</b>

<b>6</b>	<b>Battery (recharging)</b> .....	82
<b>7</b>	<b>Troubleshooting</b> .....	90
<b>8</b>	<b>Detailed information about your system</b> .....	102
<b>9</b>	<b>Index</b> .....	112

## Before you begin

Before using this product to test your blood glucose, carefully read this Owner's Booklet, and the inserts that come with the OneTouch® Verio® Test Strips and OneTouch® Verio® Control Solutions.

**NOTE:** You must charge the battery before using the OneTouch® Verio® IQ Meter. See pages 82–89 for instructions.

## IMPORTANT SAFETY INSTRUCTIONS

- This meter and lancing device are intended to be used by a single person and should not be shared, even with a family member.
- After use and exposure to blood, all parts of this kit are considered biohazardous. A used kit can potentially transmit infectious diseases even after you have performed cleaning and disinfection. For more information see:

FDA Public Health Notification: “Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens: Initial Communication” (2010)

<http://www.fda.gov/MedicalDevices/Safety/AlertsandNotices/ucm224025.htm>.

CDC Clinical Reminder: “Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens” (2010)

<http://www.cdc.gov/injectionsafety/Fingerstick-DevicesBGM.html>.

## **Intended use**

The OneTouch® Verio®IQ Blood Glucose Monitoring System is intended to be used for the quantitative measurement of glucose (sugar) in fresh capillary whole blood samples drawn from the fingertips. The system is intended to be used by a single patient and should not be used for testing multiple patients.

The OneTouch® Verio®IQ Blood Glucose Monitoring System is intended for self-testing outside the body (*in vitro* diagnostic use) by people with diabetes at home as an aid to monitor the effectiveness of diabetes control.

The OneTouch® Verio®IQ Blood Glucose Monitoring System is not to be used for the diagnosis of or screening of diabetes or for neonatal use.

Not for use on critically ill patients, patients in shock, dehydrated patients or hyperosmolar patients.

## **Test principle**

Glucose in the blood sample mixes with the enzyme FAD-GDH (refer to page 104) in the test strip and a small electric current is produced. The strength of this current changes with the amount of glucose in the blood sample. Your meter measures the current, calculates your blood glucose level, displays the blood glucose result, and stores it in its memory.

This page left blank intentionally.

# The OneTouch® Verio® IQ Blood Glucose Monitoring System

Included with your kit:

Available separately:

- A** OneTouch® Verio® IQ Meter  
(rechargeable battery included)
- B** Carrying Case
- C** OneTouch® Delica® Lancing device
- D** Sterile lancets
- E** AC Adapter
- F** Mini USB Cable
- G** OneTouch® Verio® Level 3 (Mid) Control Solution\*
- H** OneTouch® Verio® Level 4 (High) Control Solution\*
- I** OneTouch® Verio® Test Strips\*

**NOTE:**

- The OneTouch® Delica® Lancing Device uses ONLY OneTouch® Delica® Lancets.

If another type of lancing device was included, see the separate instructions for that lancing device.

\*OneTouch® Verio® Control Solutions and Test Strips are available separately. Ask for control solution where you obtain your test strips.

**⚠ WARNING**

Keep the meter and testing supplies away from young children. Small items such as the test strips, lancets, protective covers on the lancets, and control solution vial cap are choking hazards. **Do Not** ingest or swallow any items.

# 1

## Setting up your system

### Getting to know your OneTouch® Verio® IQ Blood Glucose Monitoring System

#### Meter

##### Test strip port

Insert test strip here to turn meter on

##### Display

##### Battery icon

Indicates battery power remaining

**mg/dL** is the pre-set unit of measure and cannot be changed

##### OK button

- Turns meter on/off
- Confirms menu selections

##### Back button

- Go back to the previous screen (press and release)
- Go to the Main Menu (press and hold)

##### Up and down buttons

Highlight the desired selection

##### Data/Battery charging port


- Connect to download to computer
- Connect to recharge meter battery

# Setting up your system

1

## Test strip
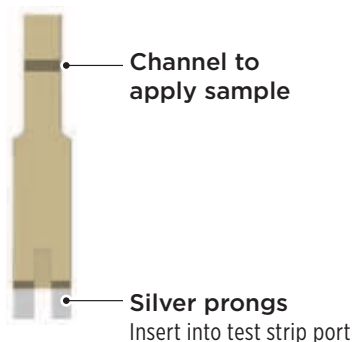
## Turning your meter on

Press and hold **OK** until the start-up screen appears. Once the start-up screen is displayed, release **OK**.

### **⚠ CAUTION**


If you see any missing pixels within the start-up screen, there may be a problem with the meter. Contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

## 1 Setting up your system

After the start-up screen, the Main Menu is displayed.

A blue bar highlights the current selection on the meter display.

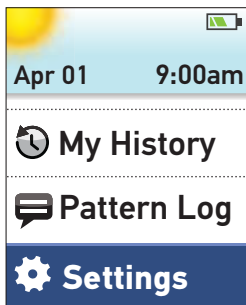**Example**  
Main Menu with  
**Settings** highlighted

## Using the display backlight

The backlight comes on automatically whenever the meter is turned on. After a few seconds of no activity, the backlight will dim. Pressing any button or inserting a test strip will turn the backlight back on.

When testing in the dark, press and hold **OK** before you insert the test strip. This will light the entire top of the meter and provides enough light to test in the dark.

## Basic settings

### Setting the meter language, time and date

Before using your meter for the first time, you should check these settings. When you turn on the meter for the first time, you will automatically start at the **Set Language** screen.

#### 1. Choose a language

In the **Set Language** screen press  or  to highlight the language you want and press .

Pressing  after making your selection confirms each setting and takes you to the next screen.

When setting up your meter you can press  to return to the previous screen to adjust a setting.

## 1 Setting up your system

**NOTE:** The time and date formats are pre-set. You cannot change these settings.

### 2. Set the time of day

In the **Set Time** screen, press ▲ or ▼ to set the hour and press OK. Repeat this step to set the minutes. Press ▲ or ▼ to set am or pm and press OK.

Set Time
12 : 45 pm
Please set with care.

### 3. Set the date

In the **Set Date** screen, press ▲ or ▼ to change the month and press OK. Repeat these steps to set the day, then the year.

Set Date
Apr 23 2011
Please set with care.

## 4. Confirm your settings

If your settings are correct press  or  to highlight **Done** and press .

If a setting is not correct, highlight the incorrect setting and press  to return to that setting screen.

Settings	
Time:	12:45pm
Date:	Apr 23 2011
Language:	English
Tool Settings	
Meter Info	
<b>Done</b>	

After confirming your settings, this screen will appear on the meter display.

To set up the High/Low Pattern tool, highlight **Set Up Now** and press .

To go to the Main Menu, highlight **Later** and press .

Settings
<b>Set up tools?</b>
Your meter is ready for testing. Set up High/Low Pattern Tool?
<b>Set Up Now</b>
Later

**NOTE:** If you select **Later**, you can always set this feature up at another time by choosing **Settings** from the Main Menu (see page 24).

# 1 Setting up your system

## Tool settings

The meter has several features that provide information that can help you manage your diabetes.

**High/Low Patterns** – Lets you know if your blood glucose levels are developing a high or low pattern.

**Tagging** – Blood glucose results can be tagged to indicate tests taken around mealtimes.

Tool Settings	
Tagging:	On
High/Low Patterns	
Alerts:	On
High/Low Limit	
Done	

## Set the High/Low Patterns

The meter can let you know whenever your blood glucose results develop a pattern below the Low limit or above the Before Meal High limit you set in your meter. See pages 66–72 for more information.

**NOTE:**

- Before Meal tags are used to determine if a blood glucose result should trigger a high pattern message. Be sure to tag your results carefully (see pages 55–59).

- Use care when adding tags to your blood glucose results. Incorrect tagging can cause the information displayed in high pattern messages to be inaccurate and/or misinterpreted.
- The High/Low Patterns feature uses your meter's time and date settings to identify patterns. Be sure the time and date are set correctly when using this feature.

## CAUTION

Be sure to talk to your health care professional about the High and Low pattern limits that are right for you. When selecting or changing your pattern limits, you should consider factors such as your lifestyle and diabetes therapy. Never make significant changes to your diabetes care plan without consulting your health care professional.

Highlight **Alerts** on the **Tool Settings** screen and press .

Tool Settings	
Tagging:	On
High/Low Patterns	
 Alerts:	On
 High/Low Limit	
Done	

## 1 Setting up your system

To turn this feature on or off, press  or  to highlight your selection and press .

Alerts	
	On
	Off
Track patterns and receive alerts.	

If the feature is turned on, highlight **High/Low Limit** on the **Tool Settings** screen and press .

### **NOTE:**

- To set high and low limits, **Alerts** must be turned on.
- Low and Before Meal High limits apply to every hour of the day.

Tool Settings	
Tagging:	On
High/Low Patterns	
 Alerts:	On
	High/Low Limit
Done	

On the **High/Low Limit** screen, highlight **Low Limit** and press **OK**.

High/Low Limit	
Low Limit:	70
High Limit:	130

The **Low Limit** is highlighted. Press **▲** or **▼** to select a value for your low limit and press **OK**.

Low Limit	
70 ▼	mg/dL
Set to track and get alerts for low glucose patterns.	

## 1 Setting up your system

Next, on the **High/Low Limit** screen, highlight **High Limit** and press **OK**.

The **High Limit** is highlighted. Press **▲** or **▼** to select a value for your high limit and press **OK**.

High Limit
Before Meal:
<b>130</b> <b>▲▼</b> mg/dL
This only applies to glucose results with Tag "Before meal".

**NOTE:** Tagging must be turned on to set a Before Meal High limit (see page 23). The meter is set at the factory with the tagging feature turned on. If you try to set the **High Limit** with the tagging feature turned off, you will be prompted to turn it on, before setting your high limit.

You will return to the **Tool Settings** screen.

Tool Settings
Tagging: On
High/Low Patterns
 Alerts: On
 High/Low Limit

## Turn the Tagging feature on or off

After testing, you can add tags to your blood glucose results. See pages 55–59 for more information.

Highlight **Tagging** on the **Tool Settings** screen and press .

Tool Settings	
<b>Tagging:</b>	<b>On</b>
High/Low Patterns	
 <b>Alerts:</b>	<b>On</b>
 <b>High/Low Limit</b>	

**NOTE:** “Before meal” tags are used to determine if a blood glucose result should trigger a high pattern message. If you want to receive high pattern messages, be sure the tagging feature is turned on and tag your results carefully (see pages 55–59).

To turn this feature on or off, press  or  to highlight your selection and press .

You will return to the **Tool Settings** screen.

Press  to return to the **Settings** screen or press and hold  to return to the Main Menu.


Tagging
<input checked="" type="checkbox"/> On
<input type="checkbox"/> Off
Turn tagging On to label glucose results Before meal or After meal.

## 1 Setting up your system

### Adjusting settings after setting up the meter for the first time


At any time, you can adjust the settings in your meter.

Highlight **Settings** on the Main Menu screen and press **OK**.

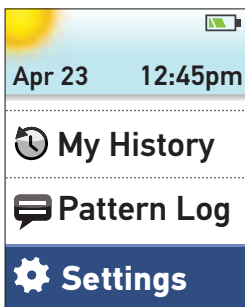Next, select **Time**, **Date**, **Language** or **Tool Settings** and press **OK**.

Then, follow the instructions on pages 15–23 to make changes to your settings.
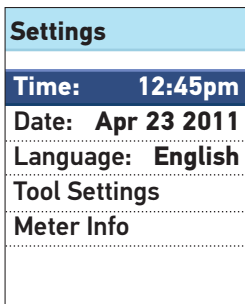
## Checking the meter serial number and software version

The meter serial number and software version are stored in your meter and you can check this information at any time.

Highlight **Meter Info** on the **Settings** screen and press **OK**.

Settings	
Time:	9:00am
Date:	Apr 01 2011
Language:	English
Tool Settings	
Meter Info	

The information is displayed.  
Press **←** or **OK** to return to the **Settings** screen.

Meter Info
Serial # <b>XXXXXXXXXX</b>
Software: 00.00
Units: mg/dL

## **1** Setting up your system

### **Turning your meter off after setup**

There are two ways to turn your meter off:

- Press and hold  for several seconds until the meter turns off.

**Or,**

- Your meter will turn off by itself if left alone for two minutes.

## Testing your blood glucose

**NOTE:** Many people find it helpful to practice testing with control solution before testing with blood for the first time. See Testing with control solution, pages 47–54.

### Preparing for a test

#### Have these things ready when you test:

OneTouch® Verio®IQ Meter

OneTouch® Verio® Test Strips

Lancing device

Sterile lancets

- Use only OneTouch® Verio® Test Strips.
- Unlike some blood glucose meters, no separate step to code your OneTouch® Verio®IQ System is required.
- Make sure your meter and test strips are about the same temperature before you test.
- **Do Not** test if there is condensation (water build-up) on your meter. Move your meter and test strips to a cool, dry spot and wait for the meter surface to dry before testing.
- Tightly close the cap on the vial immediately after use to avoid contamination and damage.
- Store unused test strips only in their original vial.
- **Do Not** open the test strip vial until you are ready to remove a test strip and perform a test. Use the test strip **immediately** after removing it from the vial.
- **Do Not** return the used test strip to the vial after performing a test.

## 2 Taking a test

- **Do Not** re-use a test strip that had blood or control solution applied to it. Test strips are for single use only.
- With clean, dry hands, you may touch the test strip anywhere on its surface. **Do Not** bend, cut or modify the test strip in any way.
- When you first open a vial of test strips, record the discard date on the label. Refer to the test strip insert or vial label for instructions on determining the discard date.

**NOTE:** Comparing your blood glucose test results taken with this meter to your results taken from a different meter is not recommended. Results may differ between meters and are not a useful measure of whether your meter is working properly. To check your meter accuracy, you should periodically compare your meter results to those obtained from a lab. Refer to page 103 for more information.

Acetaminophen, ascorbic acid (vitamin C), and uric acid (when occurring in blood at high normal or at high therapeutic concentrations) do not significantly affect results. However, abnormally high concentrations in blood may cause inaccurately low results. Tolazamide levels up to 7.5 mg/dL (3x the level found at high therapeutic levels) do not affect the blood glucose results.

**Important:** If another person assists you with testing, the meter, lancing device and caps should always be cleaned and disinfected prior to use by that person. See Care and maintenance, page 76.

### CAUTION


- The OneTouch® Verio®IQ Blood Glucose Monitoring System should not be used for patients within 24 hours of receiving a D-xylose absorption test as it may cause inaccurately high results.

### CAUTION

- **Do Not** use your test strips if your vial is damaged or left open to air. This could lead to error messages or inaccurate results. Contact Customer Service at 1 888 567-3003 (English) (Available 7 days a week, 8 a.m. – 10 p.m. Eastern Time) immediately if the test strip vial is damaged.
- If you cannot test due to a problem with your testing supplies, contact your health care professional. Failure to test could delay treatment decisions and lead to a serious medical condition.
- This meter meets the requirements for immunity to electrical interference at the frequency range and test levels specified in international standard ISO 15197:2003(E). However, testing indicates that a cell phone close to your meter during the glucose test may cause an inaccurate result or a meter Error message. **Do Not** use this meter near cell phones or similar sources of electromagnetic radiation.
- The test strip vial contains drying agents that are harmful if inhaled or swallowed and may cause skin or eye irritation.
- **Do Not** use test strips after the expiration date (printed on the vial) or the discard date, whichever comes first, or your results may be inaccurate.

## 2 Taking a test

### OneTouch® Delica® Lancing Device

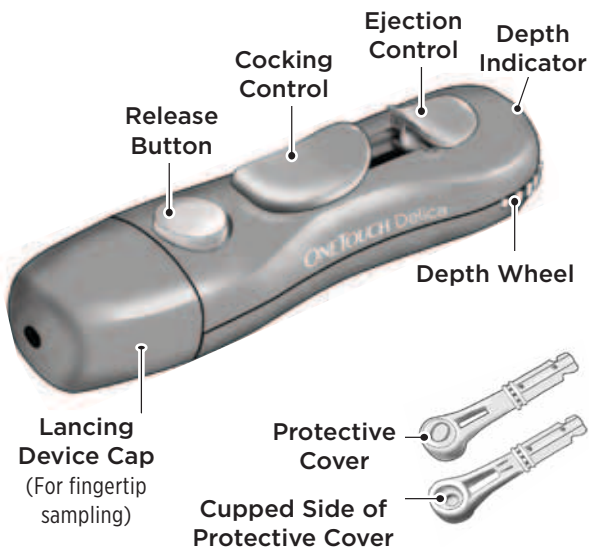**NOTE:** The OneTouch® Delica® Lancing Device uses ONLY OneTouch® Delica® Lancets.

If the lancing device shown here is different from the device included in your kit, please see the separate insert for your lancing device.

**⚠ CAUTION**

To reduce the chance of infection and disease spread by blood:

- Make sure to wash the sample site with soap and warm water, rinse and dry before sampling.
- The lancing device is intended for a single user. Never share a lancet or lancing device with anyone.
- Always use a new, sterile lancet each time you test.
- Always keep your meter and lancing device clean (see pages 76–77).

## 2 Taking a test


### Getting a blood sample from the fingertip

Choose a different puncture site each time you test. Repeated punctures in the same spot may cause soreness and calluses.

Before testing, wash your hands thoroughly with warm, soapy water. Rinse and dry completely.


#### 1. Remove the lancing device cap

Remove the cap by turning it counterclockwise and then pulling it straight off of the device.
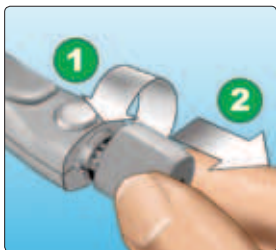
#### 2. Insert a sterile lancet into the lancing device

Align the lancet as shown here, so that the lancet fits into the lancet holder. Push the lancet into the device until it snaps into place and is fully seated in the holder.


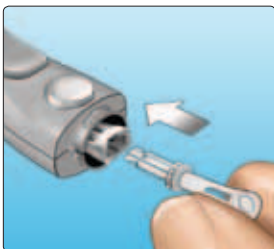
Twist the protective cover one full turn until it separates from the lancet. **Save the protective cover for lancet removal and disposal.** See pages 42–44.


### 3. Replace the lancing device cap

Place the cap back onto the device; turn clockwise to secure the cap.

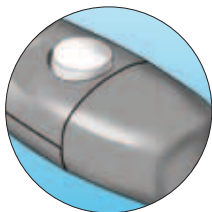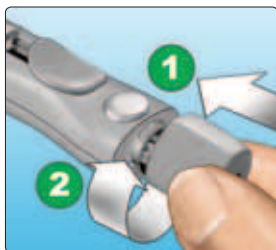
**Do Not** overtighten.
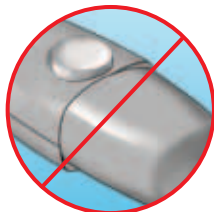
## 2 Taking a test

### 4. Adjust the depth setting

The lancing device has seven puncture depth settings, numbered 1 through 7. Smaller numbers are for a shallower puncture and the larger numbers are for a deeper puncture. Shallower punctures work for children and most adults. Deeper punctures work well for people with thick or callused skin. Turn the depth wheel to choose the setting.

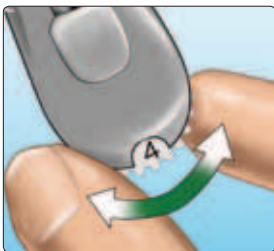**NOTE:** A shallower fingertip puncture may be less painful. Try a shallower setting first and increase the depth until you find the one deep enough to get a blood sample of the proper size.

## 5. Cock the lancing device

Slide the cocking control back until it clicks. If it does not click, it may already have been cocked when you inserted the lancet.


## 6. Insert a test strip to turn the meter on

Insert a test strip into the test strip port with the gold side of the test strip and the two silver prongs facing you.

No separate step to code the meter is required.

Silver prongs

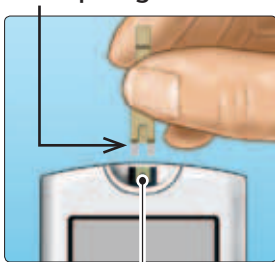
Test strip port

**NOTE:** When testing in low light conditions or in the dark, press and hold  before you insert the test strip to turn on the display backlight and test strip port light. This additional light may help you insert the test strip and complete the test.

## 2 Taking a test

When the **Apply Blood** screen appears on the display, you can apply your blood sample to either side of the test strip.

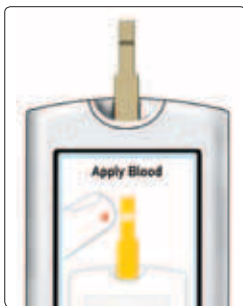
## 7. Puncture your finger

Hold the lancing device firmly against the side of your finger. Press the release button. Remove the lancing device from your finger.

## 8. Get a round drop of blood

Gently squeeze and/or massage your fingertip until a round drop of blood forms on your fingertip.

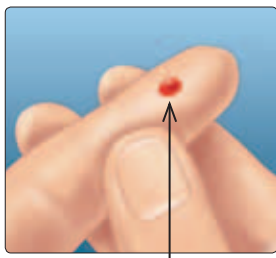Approximate size

If the blood smears or runs, **Do Not** use that sample. Dry the area and gently squeeze another drop of blood or puncture a new site.

## 2 Taking a test


### Applying blood and reading results

#### 1. Apply the sample to the test strip

You can apply blood to either side of the test strip.


Apply your sample to the opening of the channel.

Be sure to apply your sample immediately after you get a drop of blood.

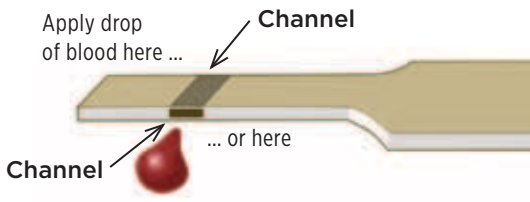
Holding the meter at a slight angle, guide the channel to the blood drop.

Channel

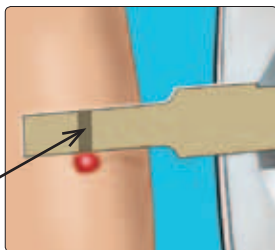When it touches your sample, the test strip wicks blood into the channel.

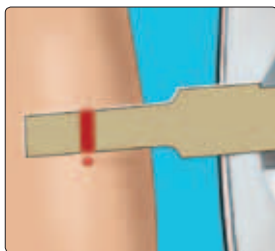
## 2 Taking a test

### 2. Wait for the channel to fill completely

The blood drop will be drawn into the narrow channel.  
The channel should fill completely.

The channel turns red and the meter will count down from 5 to 1.

Blood should **not** be applied on the top of the test strip or to the top edge of the test strip.

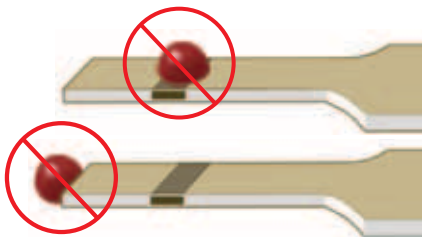- **Do Not** smear or scrape the sample with the test strip.
- **Do Not** press the test strip too firmly against the puncture site or the channel may be blocked from filling properly.
- **Do Not** apply more blood to the test strip after you have moved the drop of blood away.

- **Do Not** move the test strip in the meter during a test or you may get an error message or the meter may turn off.
- **Do Not** remove the test strip until the result is displayed or the meter will turn off.
- **Do Not** test while the battery is charging.

### 3. Read your result on the meter

Your blood glucose result appears on the display, along with the unit of measure, and the date and time of the test.

If mg/dL does not appear with the blood glucose result, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

Example

#### **⚠ CAUTION**

If the words **Control Solution** appear on the screen when testing your blood glucose, repeat the test with a new test strip. If the problem persists, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

## 2 Taking a test

### After getting a blood glucose result

Once you have your blood glucose result, you may:

- Attach a Tag to this result if the Tagging feature is turned on (see pages 55–59).

Or,

- Press and hold  to return to the Main Menu.

Or,


- Press and hold  for several seconds until the meter turns off. The meter will also automatically turn off if left alone for two minutes.

### Removing the used lancet

**NOTE:** This lancing device has an ejection feature, so you do not have to pull out the used lancet.

#### 1. Remove the lancing device cap

Remove the cap by turning it counterclockwise and then pulling it straight off of the device.
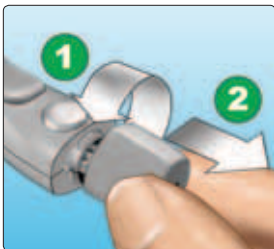
## 2. Cover the exposed lancet tip

Before removing the lancet, place the lancet protective cover on a hard surface then push the lancet tip into the cupped side of the cover.

## 3. Eject the lancet

Slide the ejection control forward until the lancet comes out of the lancing device. Return the ejection control to its back position. If the lancet fails to eject properly, cock the device again and then slide the ejection control forward until the lancet comes out.


## 2 Taking a test


### 4. Replace the lancing device cap

Place the cap back onto the device; turn clockwise to secure the cap.

It is important to use a new lancet each time you obtain a blood sample. This will help prevent infection and sore fingertips.

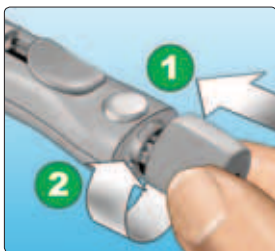**Do Not** overtighten.
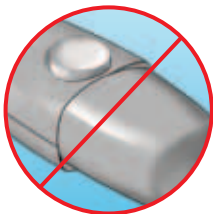
### Disposing of the used lancet and test strip

Discard the used lancet carefully after each use to avoid unintended lancet stick injuries. Used lancets and test strips may be considered biohazardous waste in your area. Be sure to follow your health care professional's recommendations or local regulations for proper disposal.

Wash hands thoroughly with soap and water after handling the meter, test strips, lancing device and cap.

## Interpreting unexpected test results

Refer to the following cautions whenever your blood glucose results are higher or lower than what you expect.

### CAUTION

#### **Low blood glucose results**

If your blood glucose result is below 70 mg/dL or is shown as EXTREME LOW GLUCOSE, it may mean hypoglycemia (low blood glucose). This may require immediate treatment according to your health care professional's recommendations. Although this result could be due to a test error, it is safer to treat first, then do another test.

#### **Dehydration and low blood glucose results**

You may get false low blood glucose results if you are severely dehydrated. If you think you are severely dehydrated, contact your health care professional immediately.

#### **High blood glucose results**

If your blood glucose result is above 180 mg/dL, it may mean hyperglycemia (high blood glucose) and you should consider re-testing. Talk to your health care professional if you are concerned about hyperglycemia.

## **2 Taking a test**

### **⚠ CAUTION**

#### **High blood glucose results**

EXTREME HIGH GLUCOSE is displayed when your blood glucose result is over 600 mg/dL. You may have severe hyperglycemia (very high blood glucose). Re-test your blood glucose level. If the result is EXTREME HIGH GLUCOSE again, this indicates a severe problem with your blood glucose control. Obtain and follow instructions from your health care professional immediately.

#### **Repeated unexpected blood glucose results**

If you continue to get unexpected results, check your system with control solution. See Control solution testing, pages 47–54.

If you are experiencing symptoms that are not consistent with your blood glucose results and you have followed all instructions in this Owner's Booklet, call your health care professional. Never ignore symptoms or make significant changes to your diabetes management program without speaking to your health care professional.

#### **Unusual red blood cell count**

A hematocrit (percentage of your blood that is red blood cells) that is either very high (above 60%) or very low (below 20%) can cause false results.

## Testing with control solution

OneTouch® Verio® Control Solution is used to check that the meter and test strips are working together properly and that the test is performing correctly. (Control solution is available separately.)

### **NOTE:**

- Use only OneTouch® Verio® Control Solution, Level 3 (Mid) or Level 4 (High), with your OneTouch® Verio® IQ Meter. Either level can be used to check your system.
- When you first open a new vial of control solution, record the discard date on the vial label. Refer to the control solution insert or vial label for instructions on determining the discard date.
- Tightly close the cap on the control solution vial immediately after use to avoid contamination or damage.

## 2 Taking a test

### CAUTION


- **Do Not** swallow or ingest control solution.
- **Do Not** apply control solution to the skin or eyes as it may cause irritation.
- **Do Not** use control solution after the expiration date (printed on the vial label) or the discard date, whichever comes first, or your results may be inaccurate.

### Do a control solution test

- Whenever you open a new vial of test strips.
- If you suspect that the meter or test strips are not working properly.
- If you have had repeated unexpected blood glucose results.
- If you drop or damage the meter.


## Performing a control solution test

1. Insert a test strip to turn the meter on

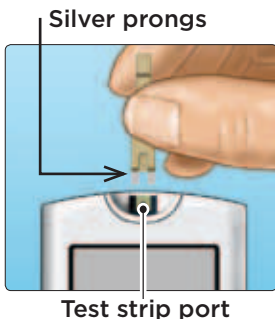Wait for the **Apply Blood** screen to appear on the display.

**NOTE:** The same **Apply Blood** screen that appears during a blood glucose test also appears during a control solution test.
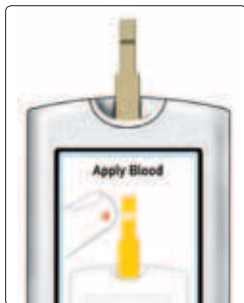
## 2 Taking a test

### 2. Prepare the control solution

Remove the vial cap and place it on a flat surface with the top of the cap pointing up.

Squeeze the vial to discard the first drop.

#### Example

Level 3 (Mid) Control Solution

Wipe both the tip of the control solution vial and the top of the cap with a clean, damp tissue or cloth.

Then, squeeze a drop into the small well on the top of the cap or onto another clean, non-absorbent surface.

### 3. Apply the control solution

Hold the meter so that the side edge of the test strip is at a slight angle to the drop of control solution.

Touch the channel on the side of the test strip to the control solution. Wait for the channel to fill completely.

## 2 Taking a test

### 4. Read your result

The meter will count down from 5 to 1. Your result is displayed along with the date, time, unit of measure, and the words **Control solution**.

The meter automatically marks the result as a control solution test.

Example

### ⚠ CAUTION

If the words **Control solution** do not appear on the screen, this result will be included in your averages and your averages will change too. Repeat the test with a new test strip. If the problem persists, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

## 5. Check if the result is in range

Each vial of test strips has both OneTouch® Verio® Level 3 (Mid) and OneTouch® Verio® Level 4 (High) Control Solution ranges printed on its label.

Compare the result displayed on the meter to either the OneTouch® Verio® Level 3 (Mid) and OneTouch® Verio® Level 4 (High) Control Solution range printed **on the test strip vial**, depending on the type of control solution you used.

Level 3 (Mid) Control Range	102–138 mg/dL
Level 4 (High) Control Range	298–403 mg/dL

### Example ranges

## Out-of-range results may be due to:

- Not following the instructions detailed on pages 49–53.
- Control solution is contaminated, expired, or past its discard date.
- Test strip or test strip vial is damaged, expired, or past its discard date.
- Meter, test strips and/or control solution were not all at the same temperature when the control solution test was performed.
- A problem with the meter.
- Dirt or contamination in the small well on the top of the control solution cap (see Step 2).

## 2 Taking a test

### 6. Cleaning

Clean the top of the control solution cap with a clean, damp tissue or cloth.

Control solution results can be seen when reviewing past results, but are not included in result averages.

#### CAUTION

- If you continue to get control solution results that fall outside the range printed on the test strip vial, **Do Not** use the meter, test strips, or control solution. Contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
- The control solution ranges printed on the test strip vial are for control solution tests only **and are not recommended ranges for your blood glucose level.**

## Attaching tags to your blood glucose results

3

Your OneTouch® Verio®IQ Meter allows you to attach a Before meal or After meal tag to your current blood glucose result.

- A Before meal blood glucose test is taken just before the start of your meal.
- An After meal blood glucose test is typically taken one to two hours after the start of your meal.

You can also change or add a tag to your most recent result (Last Result) when reviewing it in the meter memory.

We suggest you talk to your health care professional to see how meal tags may help you manage your diabetes.

**NOTE:** A Before meal tag is used to determine if a blood glucose result should trigger a high pattern message (see pages 66–72).

Use care when adding tags to your blood glucose results. Incorrect tagging can cause the information displayed in averages and high pattern messages to be inaccurate and/or misinterpreted.

### 3 Attaching tags to your blood glucose results

## Adding or changing a tag

To add tags, the Tagging feature must be turned on (see page 23).

### Adding a tag

After testing and your blood glucose result is displayed, **Before meal** and **After meal** appear on the meter display.

**Before meal**

**After meal**


## Attaching tags to your blood glucose results

3


### 1. Press or to highlight Before meal or After meal and press

If you do not want to tag this blood glucose result, simply press  to return to the result screen or press and hold  to go to the Main Menu.

**NOTE:** Once you add a tag to a blood glucose result, you may change it, but a tag cannot be removed from the result (see pages 58–59).

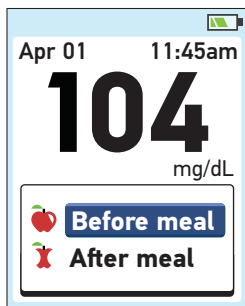Your tag will appear here →

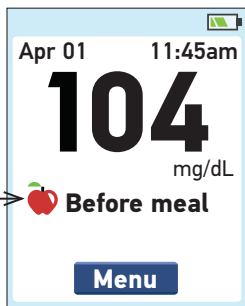
Example  
Before meal


### 3 Attaching tags to your blood glucose results

#### Adding or changing a tag on your last result

Choose **My History** from the Main Menu and press **OK** (see Reviewing past results, pages 60–61).

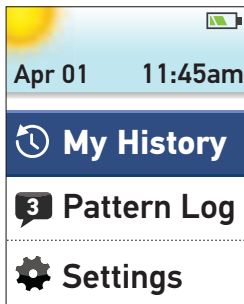Next highlight **Last Result** and press **OK**.

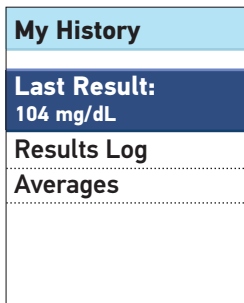
## Attaching tags to your blood glucose results


3

To add a tag highlight **Before meal** or **After meal** and press **OK**.

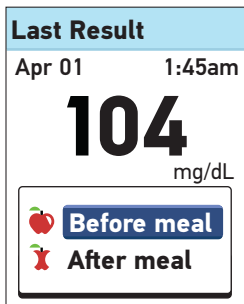To change a tag, highlight the existing tag, press **OK** and follow step 1 on page 57 to change the tag.

**NOTE:** You may only add or change a tag on your most recent (last) blood glucose result.

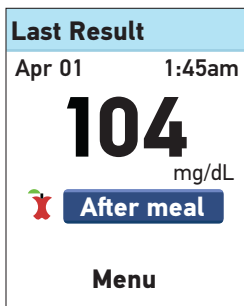
**Example**  
After meal

## 4

## Reviewing past results, averages and high/low pattern messages

Your meter stores your most recent 750 blood glucose and control solution test results and displays the blood glucose results in several ways.

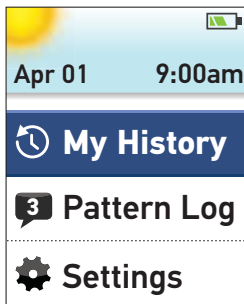If your meter is off, press and hold  until the start-up screen appears. Once the start-up screen is displayed release . (See Turning your meter on, page 13.)

If you have just completed a test, press and hold  to get to the Main Menu from the blood glucose result screen.

From the Main Menu, choose:

- **My History** to view:
  - Last Result
  - Results Log
  - Averages
- **Pattern Log** to view high or low pattern messages

Press  or  to highlight your selection and press .


# Reviewing past results, averages and high/low pattern messages

4

## My History

From the My History menu, choose:

- Last Result
- Results Log
- Averages

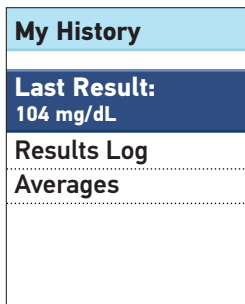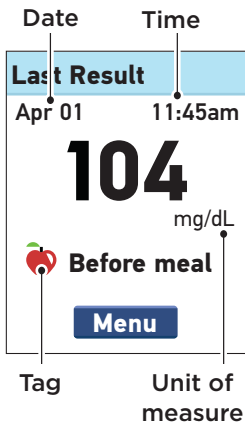
## Last Result

The meter will display your most recent result.

If you wish to add or change a tag for your Last Result, see pages 58–59.

To return to the previous screen, press .



Press  to return to the Main Menu.

## 4 Reviewing past results, averages and high/low pattern messages

### Results Log

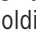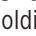The meter will display up to six results at a time, starting with the most recent.

Press  to move backward and  to move forward through your results. Pressing and holding  or  allows you to move more quickly.

Results Log		
		mg/dL
Tue, Apr 01		
9:03am		120
7:00am		110
Mon, Mar 31		
8:29pm		150
3:12pm		107

The following symbols may also appear:

**HI** if the blood glucose result was above 600 mg/dL

**LO** if the blood glucose result was below 20 mg/dL

**C** if the result is from a control solution test  
(see pages 47–54)

 if the blood glucose result is tagged Before meal

 if the blood glucose result is tagged After meal


## Reviewing past results, averages and high/low pattern messages

4

To view details of an individual result, press ▲ or ▼ to highlight the result, then press OK.

To return to the previous screen, press ↶.

Press OK to return to the Main Menu.
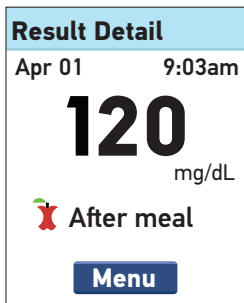
## 4 Reviewing past results, averages and high/low pattern messages

### Averages

For each of the 7, 14, 30, and 90 day periods leading up to the current date, the meter will display the number of results and the average of those results.

Averages	
	mg/dL
<b>7 Days</b> 24 Results	<b>104</b>
<b>14 Days</b> 52 Results	<b>120</b>
<b>30 Days</b> 136 Results	<b>124</b>
<b>90 Days</b> 396 Results	<b>112</b>

To view details about your 7, 14, 30, and 90 day averages, highlight the Average and press **OK**.

From the Day Avg screen, press **←** to return to the Averages screen.

Press **OK** to return to the Main Menu.

7 Day Avg
In the last 7 days, you have had 4 Low results and 5 High results tagged Before meal.
<b>Menu</b>

## Reviewing past results, averages and high/low pattern messages

4

If you do not have results in the past 7, 14, 30, and 90 day periods, the number next to Results will be zero and dashes will appear in the mg/dL column.

In result averages, an EXTREME HIGH GLUCOSE result is always counted as 600 mg/dL, and an EXTREME LOW GLUCOSE result is always counted as 20 mg/dL. (See pages 45–46 for more information on high and low blood glucose results.)

**NOTE:** The meter calculates averages based on the 7, 14, 30, and 90 day periods ending on the current date setting. If you change your date setting, your averages may change too.

Result averages provide information from past results.

**Do Not** use result averages to make immediate treatment decisions. Always consult your health care professional before making significant changes to your diabetes care plan.

### CAUTION

**Do Not** allow other people to use your meter as it may affect your averages.

## **4** Reviewing past results, averages and high/low pattern messages

### **Pattern Log (high/low pattern messages)**

High and low pattern messages appear when the meter identifies a pattern of glucose results that fall outside the High and Low limits you set on pages 18–22. Pattern messages provide information from past results and may help you to identify a need to consider making therapy or lifestyle adjustments. Always consult your health care professional before making significant changes to your diabetes care plan.

Every time you test your blood glucose, your OneTouch® Verio® IQ Meter searches for any new patterns that have developed over the past 5 days. Patterns are identified based on the time of day the test was taken. For a group of results to be considered a pattern, the time of day for each result must be within 3 hours of each other.

- A high pattern will only include results tagged “Before meal.” A high pattern message appears when your meter records 3 results that are over the High limit set in your meter.
- A low pattern message appears when your meter records 2 results below the Low limit set in your meter. No tagging is needed to trigger a low pattern.

## Reviewing past results, averages and high/low pattern messages

4

To receive high or low pattern messages with your blood glucose results, the High/Low Pattern and Tagging features must be turned on (see pages 18–23).

- Once a result is used in a pattern it will not be used again in future pattern messages.

**NOTE:** To be sure that high/low pattern messages appear whenever appropriate:

- Be sure the time and date are set correctly and update if you change time zones.
- Test your blood glucose using only this meter. Using different meters may cause you to miss patterns.
- Tag your blood glucose results carefully.
- Test when you are feeling high or low.

### CAUTION

- **Do Not** use messages to make immediate and/or significant changes to your diabetes care plan without first consulting your health care professional.
- **Do Not** wait for pattern messages to treat low or high results. Always use your current result for immediate treatment decisions.

## 4 Reviewing past results, averages and high/low pattern messages

### ⚠ CAUTION


- High and low pattern messages are based on values you set in the meter (see pages 18–22). These **Messages** are different from **EXTREME HIGH GLUCOSE** and **EXTREME LOW GLUCOSE Warnings** that appear whenever your blood glucose level is above 600 mg/dL or below 20 mg/dL.
- **Do Not** allow other people to use your meter as your patterns may be affected.

After a blood glucose test, **Low Pattern** or **High Pattern** and the message icon (💬) will appear with your test result, whenever a pattern is detected. Inaccurate results may cause pattern messages to appear.


To view a low or high pattern message after a test, highlight **View** and press (OK).

If you do not want to view the pattern message highlight **Later** and press (OK) to return to the result screen.

**NOTE:** If you choose **Later**, the flashing message icon (💬) continues to appear after each test for 48 hours or until you view your message.

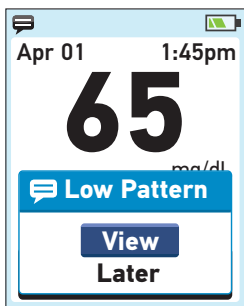Example  
Low Pattern

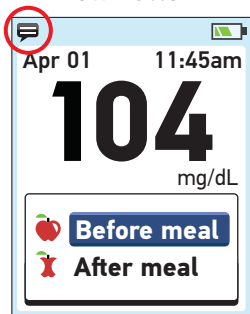
## Reviewing past results, averages and high/low pattern messages

4

### Low Pattern Message

In this example, the Low Pattern message indicates that you had at least one blood glucose result on any 2 of the past 5 consecutive days that was below the Low limit set in your meter.

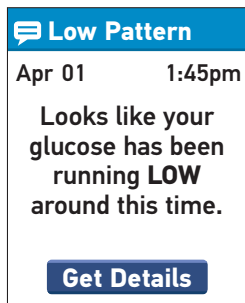Example

To view details about the individual blood glucose results that resulted in the Low Pattern message, highlight **Get Details** and press **OK**.

To return to the previous screen, press **←**.

Press **OK** to return to the Main Menu.

The screenshot shows a blue header with the text 'Low Details'. Below the header is a table with the following data:

	mg/dL
Apr 01	
1:45pm	65
Mar 30	
11:45am	60

At the bottom of the screen is a blue button with the text 'Menu'.

Example

## 4 Reviewing past results, averages and high/low pattern messages

### High Pattern Message

In this example, the High Pattern message indicates that you had at least one blood glucose result on any 3 of the past 5 consecutive days that was above the Before meal High Limit set in your meter.

 **High Pattern**

Apr 01 1:45pm

**Heads Up.**  
**Your Before meal**  
**glucose has been**  
**running HIGH**  
**around this time.**

**Get Details**

#### Example

**NOTE:** High pattern messages only include blood glucose results that are above the Before meal High limit and are tagged as Before meal results. Be sure to tag your results carefully.

To view details about the individual blood glucose results that resulted in the High Pattern message, highlight **Get Details** and press .

To return to the previous screen, press .

Press  to return to the Main Menu.

High Details		
		mg/dL
Apr 01		
1:45pm		183
Mar 31		
12:15pm		162
Mar 29		
2:03pm		190
<b>Menu</b>		


#### Example

## Reviewing past results, averages and high/low pattern messages

4

When any High or Low Pattern message has not been viewed, **Pattern Log** appears on the Main Menu along with the number of unread messages. This reminds you to view your messages.

To view the message(s), highlight **Pattern Log** and press **OK**.

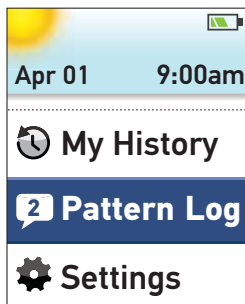On the **Pattern Log** screen, highlight the message(s) you want to view and press **OK**. After viewing a message, the message icon (🗨️) will be removed.

When you review your messages, the number that appears to the left of Pattern Log on the Main Menu is reduced by the number of messages you have read.

Pattern Log		
2 unread messages		
🗨️	Apr 01 1:45pm	Low
🗨️	Mar 29 1:00pm	High
	Mar 21 8:20am	Low
	Feb 19 6:18pm	High

## 4 Reviewing past results, averages and high/low pattern messages

To view details about the individual results that resulted in the High or Low Pattern message, highlight **Get Details** and press .

 **Low Pattern**

Apr 01 1:45pm

Looks like your glucose has been running **LOW** around this time.

**Get Details**

Example

To return to the previous screen, press .

Press  to return to the Main Menu.

**Low Details**

mg/dL

Apr 01  
1:45pm

65

Mar 30  
11:45am

 60

**Menu**

Example

### Downloading results to a computer

OneTouch® Diabetes Management Software (DMS) can store all of your records and help you spot patterns for planning meals, exercise, insulin dosing, and medication. To learn more about OneTouch® DMS and to order the software from LifeScan, visit [www.OneTouch.com](http://www.OneTouch.com) or contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. - 10 p.m. Eastern Time).

Connect only to a computer certified to UL 60950-1.

To transfer meter data, follow the instructions provided with the OneTouch® DMS to download the results from the meter.


**NOTE:** The mini USB cable used to download data to a computer is included with your kit. This cable is also used to recharge the meter battery. See pages 82–89 for more information.

## 4 Reviewing past results, averages and high/low pattern messages

Once the command to start the download is sent from the computer to the meter, the meter display will show **PC Connected** indicating that the meter is in communication mode.

**NOTE:** When the meter is connected to a computer, the meter battery will recharge. See pages 82–89 for more information.

**Do Not** insert a test strip while the meter is connected to a computer.
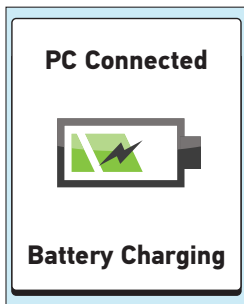
## **Reviewing past results, averages and high/low pattern messages**

**4**

This page left blank intentionally.

### Storing your system

Store your meter, test strips, control solution and other items in your carrying case. Keep in a cool, dry place between 41°F and 86°F. Keep all items away from direct sunlight and heat.

### Cleaning and disinfection

Cleaning and disinfection are different and both should be performed. Cleaning is part of your normal care and maintenance and should be performed prior to disinfection, but cleaning does not kill germs. After use and exposure to blood, all parts of this kit can potentially transmit infectious diseases. Disinfection is the only way to reduce your exposure to disease. For cleaning information, see page 77 and for disinfecting information, pages 78–81.

### Cleaning your meter, lancing device, and cap

The meter, lancing device, and cap should be cleaned whenever they are visibly dirty. For cleaning, obtain regular strength liquid dish soap and sterile gauze where you purchase household cleaning products. Prepare a mild detergent solution by stirring  $\frac{1}{2}$  teaspoon of regular strength liquid dish soap into one cup of water.

Holding the meter with the test strip port pointed down, use a sterile gauze dampened with water and mild detergent to wipe the outside of the meter and lancing device. Be sure to squeeze out any excess liquid before you wipe the meter. Wipe the inside and the outside of the cap.

Wipe dry with a clean sterile gauze.

## **5 Care and Maintenance**

### **Disinfecting your meter, lancing device, and cap**

The meter, lancing device, and cap should be disinfected periodically (at least once a week). Clean your meter, lancing device, or cap prior to disinfecting. Be sure to wipe your system completely dry as some cleaning product residue may interfere with the disinfection process.

Clorox® Germicidal Wipes\* containing 0.55% sodium hypochlorite as the active ingredient have been shown to be safe for use with the OneTouch® Verio® IQ System. Obtain them from retail websites offering disinfection products. For more information on purchase options, visit [www.OneTouchDiabetes.com/disinfection](http://www.OneTouchDiabetes.com/disinfection) or call Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. - 10 p.m. Eastern Time).

\* Other products, such as Clorox® Disinfecting Wipes, have not been tested and should not be used. Only Clorox® Germicidal Wipes should be used. Follow manufacturer's instruction for handling and storage of wipes.

Clorox® is a registered trademark of the Clorox Company.

Be sure to squeeze out any excess liquid before you wipe the meter. Hold the meter with the test strip port pointed down. Use a Clorox® Germicidal Wipe to wipe the outside of the meter, lancing device, and cap until the surface is damp.

Allow the surface of the meter, lancing device, and cap to remain damp for 1 minute.

Wipe dry with a clean, sterile gauze and allow to air dry.

Wash hands thoroughly with soap and water after handling the meter, lancing device and cap.

## 5 Care and Maintenance

- **Do Not** use alcohol or any other solvent.
- **Do Not** allow liquids, dirt, dust, blood, or control solution to enter the test strip port or the data port.
- **Do Not** squeeze the wipe or gauze into test strip port.
- **Do Not** spray cleaning solution on the meter and lancing device.
- **Do Not** immerse the meter and lancing device in any liquid.

If you have questions about cleaning or disinfecting, or if you see evidence of physical damage, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. - 10 p.m. Eastern Time).

The OneTouch® Verio®IQ System withstood cleaning and disinfection cycles well in excess of LifeScan's recommendation. See table below for more details.

	Example Usage	Tested
Weekly Cleaning	156 cleaning cycles (52 weeks x 3 years)	2879 cleaning cycles
Weekly Disinfecting	156 disinfecting cycles (52 weeks x 3 years)	412 disinfecting cycles

Though not observed in the course of testing, examples of damage to the meter may include fogged display, cracked housing or lens, illegible labels, buttons not working, or meter malfunction (such as repeated error messages). Examples of damage to the lancing device and cap may include cracking, illegible depth setting numbers, and lancing device malfunction (such as failure to load, cock or release).

## 6

# Battery (recharging)

## Battery

Your OneTouch® Verio®IQ Meter uses a rechargeable battery. When fully charged, the meter will perform blood glucose or control solution tests for about 1 to 2 weeks before recharging is needed.

**NOTE:** When the battery can no longer hold a charge, the meter needs to be replaced.

## CAUTION


The meter battery is permanent and is not replaceable.

**Do Not** try to open the meter and remove or replace the battery.

## Battery power indicator

When your meter is on, a battery icon always appears in the top right corner of the display that shows the strength of the battery.

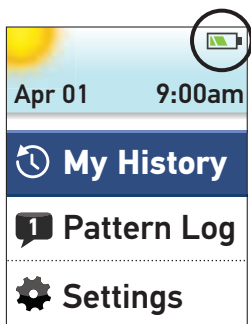**Full Power**

**30-85% Power**

**Low Power**  
(icon flashes)

**Very Low Power**  
(icon flashes)

Other screens will appear when you turn the meter on to notify you when the battery needs to be recharged. See Troubleshooting (pages 98–99) for more information.

## 6 Battery (recharging)

### Recharging the meter battery

The meter battery may be charged using one of the following options:

- Mini USB cable (computer charging)
- Mini USB cable with the AC adapter (wall charging)

Both the mini USB cable and AC adapter are included in your kit.

### CAUTION

- Only use the OneTouch® mini USB cable and AC adapter with the OneTouch® Verio® IQ Meter. If you misplace the mini USB cable or AC adapter, call Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
- **Do Not** charge meter outdoors or in a wet area.
- **Do Not** use the mini USB cable, AC adapter or meter if it is damaged, discolored, abnormally hot, or has an unusual odor. Call Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
- **Do Not** plug the AC adapter into wall outlet and leave it unattended.
- Only plug the mini USB cable included with your kit into the AC adapter.
- Verify that the wall outlet voltage matches the AC adapter voltage.
- **Do Not** allow unsupervised children to charge the meter battery.

## 1. Connect the end of the mini USB cable to the meter

Start with the meter turned off. Insert the mini USB cable into the data/battery charging port located on the bottom of your meter.

## 2. Plug the cable into the power source

Plug the other end of the cable to the USB port on your computer.

**NOTE:** Connect only to a computer certified to UL 60950-1.

**Or,**

Plug the other end of the cable into the AC adapter. Then, plug the AC adapter into the wall outlet.

## 6 Battery (recharging)

### WARNING


To avoid a possible shock, **Do Not** insert a test strip while the battery is charging.

### **NOTE:**


- Using the mini USB cable or AC adapter fully charges the batteries in about 2 hours.
- If using the USB port on your computer to charge the battery, be sure the computer is turned on and not in standby mode. If the meter does not charge, try using another USB port on your computer.
- To optimize battery life, it is best to recharge the battery when the **Low Battery** screen appears (see page 98).

### 3. Charge the battery

When the meter is connected to a computer or when using the AC adapter, **Battery Charging** appears on the display to indicate that the battery is charging.

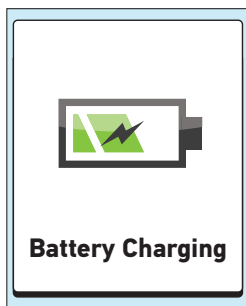When the meter is connected to a computer and you have installed OneTouch®Diabetes Management Software (see pages 73–74), **PC Connected** also appears on the meter display.

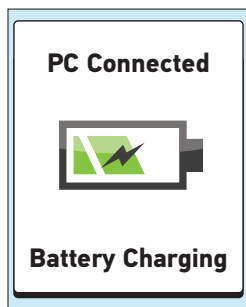**NOTE:**



- Using either charging option, **Battery Charging** must appear or the meter battery is not charging.
- After 20 seconds the **Battery Charging** screen will dim and the meter will continue to charge.

## 6 Battery (recharging)

### 4. Complete the charge


When the charge is completed (about 2 hours when charging an empty battery) the **Charge Complete** screen appears. Remove the mini USB cable from the meter before turning the meter on or before you take a test.

**NOTE:** If the meter is not fully charged after 2 hours, check that the cable is plugged in correctly. If this does not solve the problem, call Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
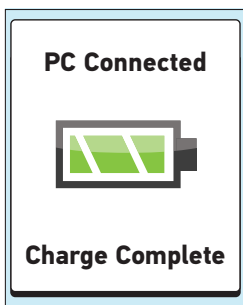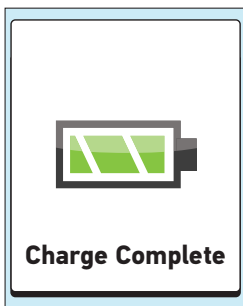
## Rapid charging

If the battery power is too low to take a blood glucose test, you can connect the meter (USB or AC adapter) for a rapid charge. **Rapid Battery Charging** appears on the display for approximately 60 seconds to indicate the meter is in rapid charge mode. After rapid charge, the battery will continue to charge if left connected.

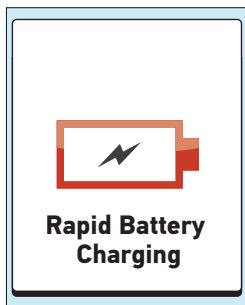After the rapid charge, be sure to disconnect the meter from the wall outlet or computer **before testing your blood glucose**. After testing, reconnect your meter and complete the charge.

**NOTE:** Even when battery power is extremely low, it will not affect results stored in the meter. However you may need to reset the date and time settings. See pages 15–23.


## Meter disposal

Dispose of the meter and internal lithium polymer battery according to your local/country laws and regulations.

## Troubleshooting

The OneTouch® Verio®IQ Meter displays messages when there are problems with the test strip, with the meter, or when your glucose levels are above 600 mg/dL or below 20 mg/dL. Improper use may cause an inaccurate result without producing an error message.

**NOTE:** If the meter is on but does not operate (locks-up), call Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
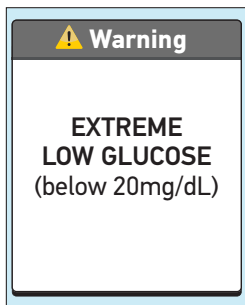
### What it means

You may have a very low blood glucose level (severe hypoglycemia), below 20 mg/dL.

### What to do

**This may require immediate treatment.** Although this message could be due to a test error, it is safer to treat first and then do another test. Always treat according to your health care professional's recommendations.

## Warning

**EXTREME  
HIGH GLUCOSE**  
(above 600mg/dL)

### What it means

You may have a very high blood glucose level (severe hyperglycemia), over 600 mg/dL.

### What to do

**Re-test your blood glucose level.** If the result is EXTREME HIGH GLUCOSE again, obtain and follow instructions from your health care professional right away.

## 7 Troubleshooting

### Warning

**Temperature  
too high. Out of  
operating range.  
See Owner's  
Booklet.**

### What it means

Meter is too hot (above 111°F) to work correctly.

### What to do

Move the meter and test strips to a cooler area. Insert a new test strip when the meter and test strips are within the operating range (43–111°F). If you do not get another Temperature too high message, you can proceed with testing.

**Warning**

**Temperature  
too low. Out of  
operating range.  
See Owner's  
Booklet.**

**What it means**

Meter is too cold (below 43°F) to work correctly.

**What to do**

Move the meter and test strips to a warmer area. Insert a new test strip when the meter and test strips are within the operating range (43–111°F). If you do not get another Temperature too low message, you can proceed with testing.

## 7 Troubleshooting

### ERROR 1

**Meter problem.  
Contact Customer  
Service.**

### What it means

There is a problem with the meter.

### What to do

**Do Not** use the meter.  
Contact Customer Service at  
1 888 567-3003 (English)  
(available 7 days a week,  
8 a.m. – 10 p.m. Eastern Time).

### ERROR 2


**Meter or strip  
problem. Retest  
with a new strip.**

### What it means

Error message could be caused  
either by a used test strip or a  
problem with the meter.

### What to do

Repeat the test with a new  
test strip; see blood  
application (pages 38–41)  
or control solution testing  
(pages 47–54). If this message  
continues to appear, contact  
Customer Service at  
1 888 567-3003 (English)  
(available 7 days a week,  
8 a.m. – 10 p.m. Eastern Time).
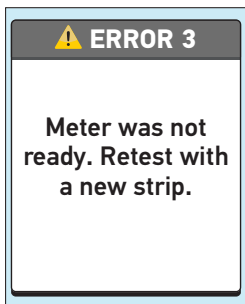
### What it means

The sample was applied before the meter was ready.

### What to do

Repeat the test with a new test strip. Apply a blood or control solution sample only after **Apply Blood** appears on the display. If this message continues to appear, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

### ERROR 4

**Strip fill problem.  
Retest with a  
new strip.**


### What it means

*One of the following may apply:*

- Not enough blood or control solution was applied or more was added after the meter began to count down.
- The test strip may have been damaged or moved during testing.
- The sample was improperly applied.
- There may be a problem with the meter.

### What to do

Repeat the test with a new test strip; see blood application (pages 38–41) or control solution testing (pages 47–54). If the error message appears again, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).
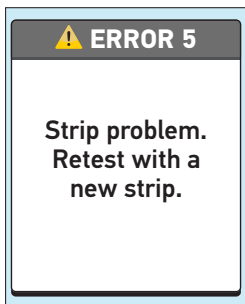
### What it means

The meter has detected a problem with the test strip. Possible cause is test strip damage.

### What to do

Repeat the test with a new test strip; see blood application (pages 38–41) or control solution testing (pages 47–54). If the error message appears again, contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time).

## 7 Troubleshooting

**Low Battery**  
Charge soon.

**OK**

### What it means

Battery power is low but there is still enough power to perform a test.

### What to do

Press **OK** to continue but recharge the battery as soon as possible.

**Very Low Battery**  
Charge now.

**OK**

### What it means

Battery power is extremely low but there is still enough power to perform a test.

### What to do

You may press **OK** to continue but you should recharge the battery immediately.


## What it means

There is not enough battery power to perform a test.

## What to do

Recharge the battery now. If you need to test your blood glucose right away, use the rapid charging option (see page 89 for more information on rapid charging).
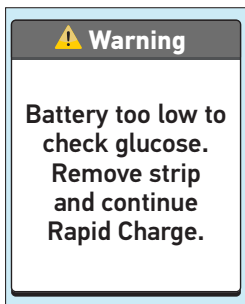
## What it means

There is not enough battery power to perform a test.

## What to do

Remove the test strip and recharge the battery. See page 89 for information on rapid charging.

## 7 Troubleshooting

### Warning

**Result not available.**

### What it means

*One of the following may apply:*

- No result in memory, such as the first time use of the meter or after downloading all data to a computer.
- Your meter was unable to recall this result. This result will not be included in result averages.

### What to do

Contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time) to report this occurrence, **unless** this is your first use of the meter. You can still perform a blood glucose test and get an accurate result.

Results Log		
		mg/dL
Tue, Apr 01		
7:03am		180
12:00pm		---
6:57am		142
Mon, Mar 29		
12:00pm		110

## What it means

Your meter was unable to recall this result. This result will not be included in result averages.

## What to do

Contact Customer Service at 1 888 567-3003 (English) (available 7 days a week, 8 a.m. – 10 p.m. Eastern Time) to report this occurrence. You can still perform a blood glucose test and get an accurate result.

## Detailed information about your system

### Comparing meter results to laboratory results

Results obtained from the OneTouch® Verio®IQ Meter and laboratory tests are reported in plasma-equivalent units. However, your meter result may differ from your lab result due to normal variation. A result from your OneTouch® Verio®IQ Meter is considered accurate when it is within  $\pm 20\%$  of the lab result.

Meter results can be affected by factors that do not affect lab results in the same way, which may cause a difference of more than  $\pm 20\%$ . Specific factors that may cause your meter result to vary from your lab result by more than  $\pm 20\%$  include:

- You have eaten recently. This can cause a result from fingertip testing to be up to 70 mg/dL higher than a lab test using blood drawn from a vein.<sup>1</sup>
- Your hematocrit is above 60% or below 20%.
- You are severely dehydrated.
- For additional information, refer to the OneTouch® Verio® Test Strip Insert.

1. Sacks, D.B.: "Carbohydrates." Burtis, C.A., and Ashwood E.R. (ed.), *Tietz Textbook of Clinical Chemistry*, Philadelphia: W.B. Saunders Company (1994), 959.

## Detailed information about your system

8

### **Guidelines for obtaining accurate meter to lab comparisons:**

#### **Before going to the lab:**

- Perform a control solution test to make sure your meter is working properly.
- **Do Not** eat for at least 8 hours before you test your blood.
- Take your meter and testing supplies with you to the lab.

### **Testing with your OneTouch® Verio®IQ Meter at the lab:**

- Test within 15 minutes of the lab test.
- Use only a fresh, capillary blood sample from your fingertip.
- Follow all instructions in this Owner's Booklet for performing a blood glucose test.

### **Comparing your meter results to those taken from another meter**

Comparing your blood glucose test results taken with this meter to your results taken from a different meter is not recommended. Results may differ between meters and are not a useful measure of whether your meter is working properly.

## **8** Detailed information about your system

### **Technical Specifications**

<b>Reported result range</b>	20–600 mg/dL
<b>Calibration</b>	Plasma-equivalent
<b>Sample</b>	Fresh capillary whole blood
<b>Sample volume</b>	0.4 µL
<b>Test time</b>	5 seconds
<b>Assay method</b>	FAD-GDH (flavin adenine dinucleotide dependent glucose dehydrogenase)
<b>Power source</b>	Rechargeable 3.7 Volt lithium polymer battery
<b>Expected meter battery life</b>	1 to 2 weeks in between charges
<b>Unit of measure</b>	mg/dL
<b>Memory</b>	750 test results

Technical Specifications

Automatic shutoff	Two minutes after last action
Size	3.46 x 1.85 x 0.47 inches
Weight	Approximately 1.66 ounces
Operating ranges	Temperature: 43–111°F  Relative humidity: non-condensing 10–90%  Altitude: up to 10,000 feet  Hematocrit: 20–60%
Battery type	Rechargeable, non-serviceable, 150mAh, 3.7 Volt DC nominal, lithium polymer battery (5V input charge voltage)
Charging current	100mA,  Direct current
Charging temperature	32–122°F

## System Accuracy

Diabetes experts have suggested that glucose meters should agree within 15 mg/dL of a laboratory method when the glucose concentration is lower than 75 mg/dL, and within 20% of a laboratory method when the glucose concentration is 75 mg/dL or higher. Samples from 100 patients were tested using both the OneTouch® Verio®IQ System and the YSI 2300 Glucose Analyzer laboratory instrument.

### System Accuracy Results for Glucose Concentrations <75 mg/dL

Percent (and number) of meter results that match the laboratory test

Within ±5 mg/dL	Within ±10 mg/dL	Within ±15 mg/dL
88.2% (45/51)	100% (51/51)	100% (51/51)

**System Accuracy Results for  
Glucose Concentrations  $\geq 75$  mg/dL**

Percent (and number) of meter results that match the laboratory test

Within $\pm 5\%$	Within $\pm 10\%$	Within $\pm 15\%$	Within $\pm 20\%$
71.1% (177/249)	94.8% (236/249)	98.0% (244/249)	100% (249/249)

**System Accuracy Results across the  
entire Glucose Range**

Percent (and number) of meter results that match the laboratory test

Within $\pm 15$ mg/dL or $\pm 20\%$
100% (300/300)

Therefore, 100% of the total results obtained with the OneTouch® Verio®IQ System achieved the goal suggested by the diabetes experts.

8

Detailed information about your system

### Regression Statistics

Samples were tested in duplicate on each of three test strip lots. Results indicate that the OneTouch® Verio®IQ System compares well with a laboratory method.

# of Subjects	# of Tests	Slope	Intercept (mg/dL)
100	300	0.95	4.81

95% CI Slope	95% CI Intercept (mg/dL)	Std. Error (S <sub>y.x</sub> ) (mg/dL)	R <sup>2</sup>
0.94 to 0.96	2.76 to 6.87	9.73	0.99

Precision

Within Run Precision (300 Venous Blood Tests)

Target Glucose (mg/dL)	Mean Glucose (mg/dL)	Standard Deviation (mg/dL)	Coefficient of Variation (%)
40	42.39	1.21	2.86
100	97.43	1.84	1.89
130	130.24	2.49	1.91
200	196.15	3.96	2.02
350	330.34	7.39	2.24

Results show that the greatest variability observed between test strips when tested with blood is 2.86% or less

Total Precision (600 Control Solution Tests)

Glucose Level Ranges (mg/dL)	Mean Glucose (mg/dL)	Standard Deviation (mg/dL)	Coefficient of Variation (%)
Low (38-62)	50.31	1.77	3.52
Normal (102-138)	117.14	2.70	2.31
High (298-403)	342.34	8.68	2.53

8

Detailed information about your system

### Fingertip Results for Glucose Concentrations <75 mg/dL

Percent (and number) of meter results that match the laboratory test

Within ±5 mg/dL	Within ±10 mg/dL	Within ±15 mg/dL
58.1% (18 of 31)	100% (31 of 31)	100% (31 of 31)

### Fingertip Results for Glucose Concentrations ≥75 mg/dL

Percent (and number) of meter results that match the laboratory test

Within ±5%	Within ±10%	Within ±15%	Within ±20%
64.9% (159 of 245)	93.5% (229 of 245)	99.2% (243 of 245)	99.6% (244 of 245)

### Guarantee

LifeScan guarantees that the OneTouch® Verio® IQ Meter will be free of defects in material and workmanship for three years, valid from the date of purchase. The guarantee extends only to the original purchaser and is not transferable.

### Electrical and safety standards

This meter complies with CISPR 11: 2009, Class B (Radiated Only). Emissions of the energy used are low and not likely to cause interference in nearby electronic equipment. The meter has been tested for immunity to Level 3 electrostatic discharge as specified in IEC 61000-4-2. This meter has been tested for immunity to radio frequency interference over the frequency range 80MHz to 2.5GHz at 3V/m as specified in IEC 61000-4-3.

**Do Not** use the equipment where aerosol sprays are being used, or when oxygen is being administered.

AC adapter.....	11, 84, 85
Averages.....	64
Backlight.....	14, 35
Basic settings.....	15
Battery (charging).....	82
Battery Empty message.....	99
Buttons on meter.....	12
Charging port.....	12, 85
Cleaning your meter, lancing device, and cap.....	76, 77
Comparing meter results to laboratory results.....	102
Control solution, discard and expiration dates.....	47, 48
Control solution, testing.....	47
Data port.....	12
Date setting.....	16
Dehydration.....	45
Disinfecting your meter, lancing device, and cap.....	76, 78
Display check.....	13
Disposal, lancets and test strips.....	44
Downloading results to a computer.....	73
Error messages.....	90
EXTREME HIGH GLUCOSE message.....	46, 91
EXTREME LOW GLUCOSE message.....	45, 90
Fingertip testing procedure.....	27
Guarantee.....	111
High/Low pattern messages.....	18, 66
Hyperglycemia.....	46, 91

Hypoglycemia.....	45, 90
Infection, reduce the chance.....	31
Intended use.....	7
Kit components.....	11
Lancing device.....	11, 30
Language, setting.....	15
Low Battery icon.....	3, 83
Low Battery message.....	98
Main Menu screen.....	14, 60
mg/dL.....	12, 41, 104
Mini USB cable.....	11, 73, 84, 85
Pattern messages.....	18, 66
Plasma calibration.....	102, 104
Results, reviewing past.....	60
Serial number.....	25
Settings.....	15
Software, OneTouch® Diabetes Management.....	73
Start-up screen.....	13, 60
Storing your system.....	76
Symbols.....	3, 62
Tags feature.....	18, 23, 55
Technical specifications.....	104
Temperature.....	27, 76, 92, 93, 105
Test strip, applying drop of blood.....	38
Test strip, discard and expiration dates.....	28, 29
Test strip, silver prongs.....	13, 35, 49

## 9 Index

Time setting.....	16
Tool settings.....	18
Turning meter off.....	26, 42
Turning meter on.....	13, 35, 49, 60
Unexpected blood glucose results.....	45
Unit of measure.....	12, 41, 104
Unusual red blood cell count (hematocrit).....	46
Very low battery message.....	98

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© 2011 Lifescan, Inc.